

Keeping Baby Warm Just Got Easier...




The NEW Solution to  
Active Infant Warming

Underbody Conductive Warming with Flexible  
Temperature Settings to Minimize Risk of Infant Cold Stress



# Multi-Layer - Multi-Use - Multi-Benefit

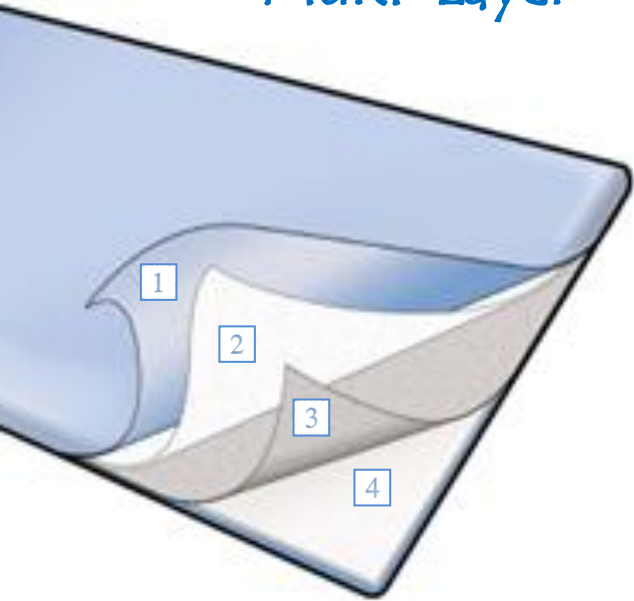

## 1 Soft, yet Durable Cover

Baby sensitive, Latex-Free. 100% reusable with welded seams for long life. Easy to clean and maintain.

## 2 Ultra Soothing Layer

Enhanced softness and comfort. Promotes even transfer of warmth without hot / cold spots.

## 3 Patented Carbon Polymer Thermal Fabric

Delivers continuous, underbody warming at specific, selected temperature setting.

## 4 Pressure Relief "Memory Foam"

Individualized for maximum support and comfort during this fragile stage of development.

Koala  $\Delta T$  Controller is user-friendly with hi-visibility control buttons for ease of use. Provides accurate transfer of information to /from mattress. Safety features include thermal cut-off to prevent overheating. Five language settings for multi-lingual capability.

Offers three temperature ranges with adjustable settings in 1.0 °F / 0.5°C. Temperature ranges offered:

91 - 99°F / 33 - 37°C  
82 - 95°F / 28 - 35°C  
90 - 102°F / 32 - 39°C

## Warming Mattress Features

- ▶ Soft, yet durable mattress design with patented conductive thermal technology.
- ▶ Delivers active underbody warmth providing a cost efficient alternative to keep baby warm.
- ▶ 100% reusable with no disposable parts. Available in a variety of sizes to fit most standard size cribs/bassinets.





## Controller Features

- ▶ Compact, energy efficient, lightweight design with hi-visibility control panel.
- ▶ Provides continuous warming with minimal warm-up time (3-5 minute).
- ▶ Temperature selection available in °F or °C.
- ▶ Adjustable pole mount for ease of use.



# KOALA<sup>ΔT</sup> Infant Warmer



PRODUCT DESCRIPTION	MODEL #	DIMENSIONS	WEIGHT
 <b>Standard Warming Mattress</b>	28 - KIW - SM1	24" x 13" 61 x 34 cm	1.3 lbs 0.6kg
 <b>Large Warming Mattress</b>	28 - KIW - LM2	29" x 23" 73 x 58 cm	2.4 lbs 1.1kg
 <b>Warming Blanket</b>	28 - KIW - BL3	18" x 20" 45 x 50 cm	0.9 lbs 0.4kg
 <b>Controller</b>	28 - KIW - CON	6.2" x 9.5" x 9.1" 16 x 24 x 23 cm	9.0 lbs 4.1kg

**NOVAMEDUSA**

4 Westchester Plaza  
 Elmsford, NY 10523  
 800-425-3535  
[www.NovaMed-USA.com](http://www.NovaMed-USA.com)